
# LAFE
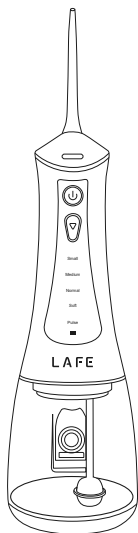
## CleanPRO

**PL**

Irygator

**EN**

Water flosser

**CS**

Irigátor

**SK**

Irygátor

**HU**

Öntözőgép

**BG**

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Model: LAFFLO47480, LAFFLO47481

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Akumulator: 2000 mAh

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Please carefully read and retain the following information for future use to ensure satisfactory and safe operation of the product.

**Safety Information:**

- This device is intended for home use only.
- Remove all packaging elements before first use.
- Unwind and straighten the power cord before use.
- Ensure the power cord has no visible damage.
- Make sure the power cord is not stretched over an open flame or other heat sources that could damage its insulation.
- The cable must be positioned to avoid the risk of tripping or entanglement.
- This product is not a toy. Ensure children are supervised and do not play with the product.
- This device is not intended for use by individuals (including children) with limited physical, sensory, or mental capabilities or lack of experience and knowledge unless supervised or instructed on the use of the device by a person responsible for their safety.
- Do not attempt to repair, modify, or disassemble any part of the device on your own.
- Do not use the device while bathing or showering.

**WARNINGS:**

- The device is not a toy.
- Cleaning and maintenance of the device should not be performed by children without adult supervision.
- The device should not be disassembled except when it is to be disposed of. The battery should be removed from the device before disposal.
- The battery should be disposed of in accordance with local regulations.
- The device should be charged only with the supplied cable.
- Do not fill the water tank with water hotter than 40°C.
- Do not strike or drop the device; discontinue use if damaged.
- Do not operate the device without an attached tip.
- Do not immerse the device in water or any other liquid.
- Keep the device tips out of reach of children to avoid the risk of choking.
- Do not pour mouthwash liquids into the water tank as they may damage the device.

- If you experience discomfort during or after using the device, stop using it and consult a doctor.
- Do not press the tips too hard against your teeth and gums.
- Use only the tips provided with the device.
- The device is intended for cleaning interdental spaces, gums, and tongue. Do not use it for any other purpose.
- Use gum pocket tips only on gentle modes. High pressure can damage the tip.

### Parts Description:

#### Fig. 1

1. Tip release button
2. ON/OFF button
3. Mode switch button
4. Charging status indicator (green - charging complete, red - connect to charge)
5. Water tank
6. USB-C port
7. Water tank lid

### Tip Types:

#### Fig. 2

- 1,3. Tips for orthodontic appliances, crowns, and bridges.
2. Gum pocket tip.
4. Tongue tip.
5. Standard tip (2 pieces included).

### USAGE:

1. Attach the selected tip to the device body.
2. Open the water tank lid and fill the device with water not exceeding 40°C.
3. Press the mode switch button to select the operating mode of the irrigator:  
Small, Medium, Normal: Modes for daily oral hygiene.  
Soft: Mode for those with sensitive gums and for gum pocket tips.  
Pulse: Mode for simultaneous cleaning and gum massage.

4. Hold the device upright with the tip aimed at your teeth and gums during use.
5. Direct the water stream at a 90° angle to the teeth and gums.
6. Slowly move the irrigator tip along the gum line and interdental spaces.
7. Pay special attention to cleaning orthodontic brackets, crowns, or bridges.
8. After irrigation, turn off the device, empty the remaining water, and dry the irrigator.

### Package Contents

Irrigator, USB-C charging cable, tips (6 pcs).

### Product Disposal

1. Paper, cardboard, and similar packaging materials should be recycled or disposed of in waste containers intended for paper.
2. Plastic elements such as styrofoam, foil, granulated fillers, and similar materials should be disposed of in containers intended for plastics.

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The symbol of a crossed-out wheeled bin placed on electronic or electrical equipment, its packaging or accompanying documents means that the product may not be thrown out together with other waste. Used equipment may contain substances with toxic and carcinogenic properties, hazardous to human health and life, and poisoning the soil and groundwater. It is the user's responsibility to hand over the used equipment to a designated collection point for its proper processing. For more information on recycling of electronic and electrical equipment, please contact your local authorities, waste disposal services and the place where you purchased this product.

The household plays an important role in contributing to the reuse and recovery of secondary raw materials, including recycling, waste equipment. At this stage, attitudes are formed that affect the preservation of the common good, which is a clean natural environment.

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